

Preparation zone



1
Clean hands
with alcohol hand rub
or soap & water



2
Clean trolley
according to
local policy



3
Gather dressing
pack & equipment,
place on bottom shelf

Patients zone




4
Apply apron
& clean hands



5
Open Critical
Aseptic Field &
equipment using NTT



6
Clean hands & apply
sterile gloves using
closed gloved technique

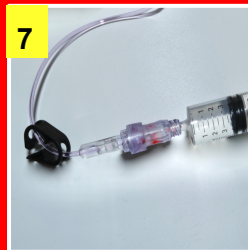

7
Prime EZ-Huber
device with 0.9%
NaCl using NTT




8
Disinfect skin
for 30 sec - 2% chlorhexidine
/70% alcohol applicator & a cross
hatch technique. Allow to dry




9
Stabilise Port-A-
Cath with non-
dominant hand



10
Insert EZ-Huber
device aiming for centre of
Port-A-Cath at a 90-degree
angle

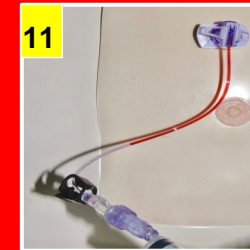

11
Confirm correct
placement with
'flashback' of blood
in extension tubing

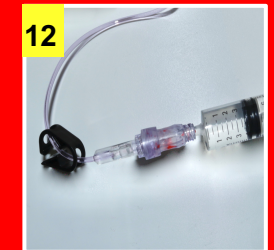

12
Flush EZ-Huber
as per local policy



13
Apply dressing
as per local policy



14
Dispose of sharps
& equipment



15
Dispose of gloves
and then apron



16
Clean hands
Immediately after
glove removal

Decontamination zone



17
Clean trolley
according to
local policy



18
Clean hands
with alcohol hand rub
or soap & water